

Easy & Safe! Sinus solution proposed by Advisory Committee member of MegaGen Implant

Advanced Sinus Solution

- Improved perforation stability + instant bone grafting
- Simple & safe hydraulic lift with 1mm step-by-step stoppers
- Precise elevation, even with minimal window opening
- Diamond drill & bone carrier exclusively for lateral approach

Simple surgery concept

MICA Kit


- Crestal Approach Sinus kit, Advanced!



Multifunctional,
Hybrid Flute


Advanced CASA Drill

Improved perforation stability + instant bone grafting

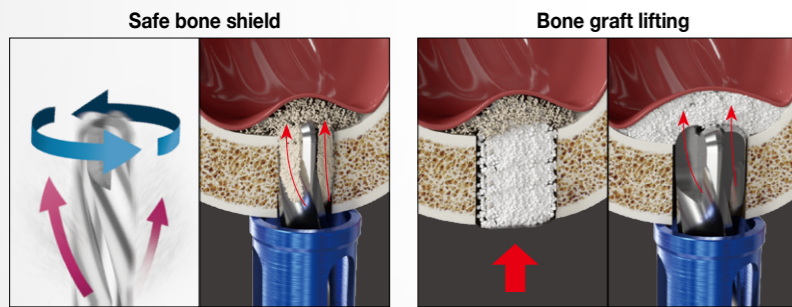

Safe bone shield
When drilling, bone chips are expelled towards the sinus floor, significantly reducing the risk of membrane perforation.

Bone graft lifting
With a low speed, bone graft lifting can be performed without any additional instruments.

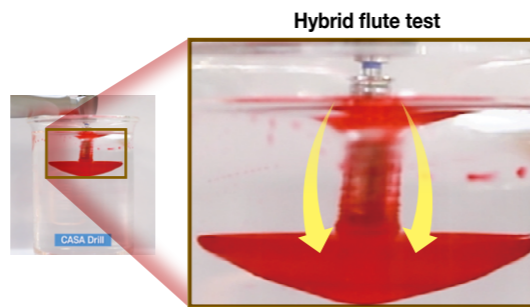

Hybrid flute test

CASA Drill

Minimalized window concept

MILA Kit

- Lateral Approach Sinus kit, Advanced!



Advanced Sinus Elevator

Precise control, even with minimal window space

The ergonomic design of the Master Grip is optimized for a lateral approach to enable meticulous and precise finger control even with a minimized window.




Depth Gauge

Elevator 1

Elevator 2

Elevator 3


Tips for step-by-step use in various environments to enable safe membrane lifting.

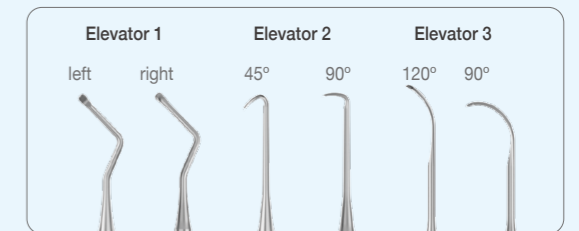

Advanced Aqua Lifter (coming soon)

Simple & safe hydraulic lift

Since a hydraulic sinus lift procedure must be completely sterile, MegaGen's single-use Aqua Lifter completely eliminates any source of instrument contamination, plus it connects directly to the syringe (no silicon tube connection), making it easy and convenient to use.

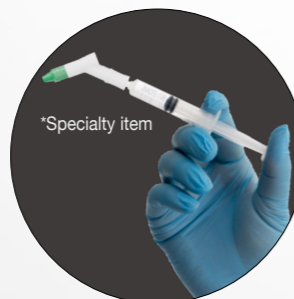

*Specialty item

Advanced Drill Stopper



2mm

3mm

4mm

5mm

6mm

7mm

8mm

9mm

10mm

11mm


12mm

Delicate depth adjustment in 1mm increments

11 stoppers with safe 1mm step-by-step increments from 2 - 12mm to accommodate a wide range of cases.

Advanced Diamond Drill

Safe & easy for opening lateral wall

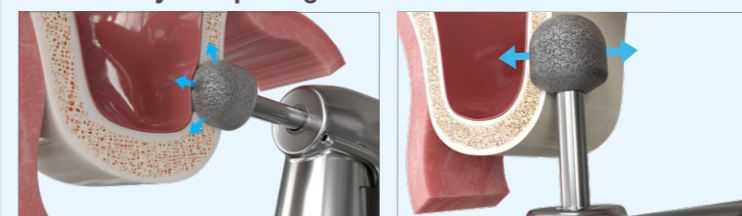

Ø6.5 diamond drill for easy and safe formation of a minimal window.

Designed to facilitate lateral cutting, making it easy to use, even in hard-to-reach situations.

Easy Bone Carrier

Instant supply of large quantity of bone graft material



MAX 0.5cc

Specially designed for bone grafting with a minimized window to enable the instant transfer of a large quantity of bone graft material.

Precise Control,
Master Grip



Simple Surgery Concept

MICA Kit


Ref.C

MKCA3000S

· Crestal Approach Sinus kit, Advanced!




Components

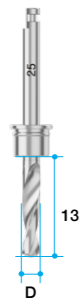

Initial Drill

Diameter	Ref.C
Ø2.5	SCID25

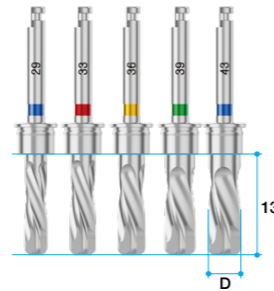


CASA Drill

Diameter	Ref.C
Ø2.9	SCAD29
Ø3.3	SCAD33
Ø3.6	SCAD36
Ø3.9	SCAD39
Ø4.3	SCAD43




Drill Stopper

Length	Ref.C
2mm	SCDS02
3mm	SCDS03
4mm	SCDS04
5mm	SCDS05
6mm	SCDS06
7mm	SCDS07
8mm	SCDS08
9mm	SCDS09
10mm	SCDS10
11mm	SCDS11
12mm	SCDS12

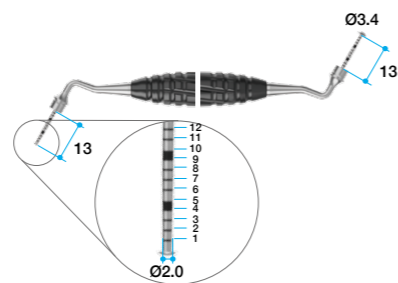

Depth Gauge

Diameter	Ref.C
Ø2.0 / Ø3.4	SCDG2034

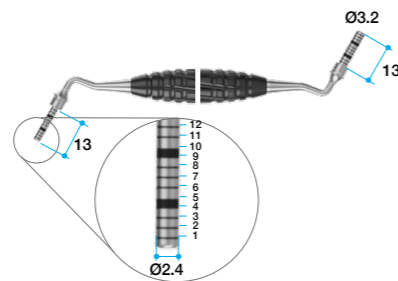

Bone Condenser

Diameter	Ref.C
Ø2.4 / Ø3.2	SCBC2432



Aqua Lifter (coming soon)

Diameter	Ref.C
Ø6.2	*SCAL62P

(*) Separate sales item.


Minimalized Window Concept

MILA Kit


Ref.C

MKLA3000M

· Lateral Approach Sinus kit, Advanced!

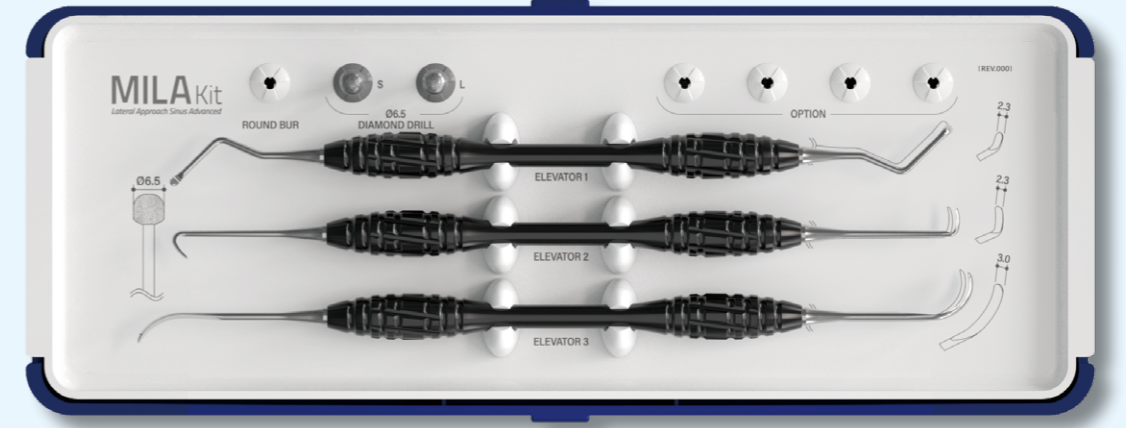


Components

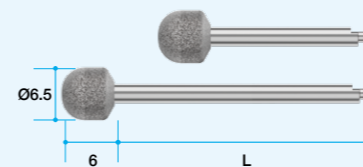

Diamond Drill

Length	Ref.C
20mm	SLDD65S
35mm	SLDD65L

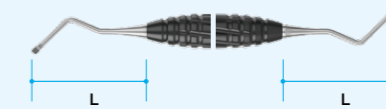


Elevator 1

Length	Ref.C
39mm	SLE100




Elevator 2

Length 1	Length 2	Ref.C
39mm	41mm	SLE200

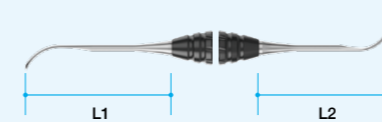


Elevator 3

Length 1	Length 2	Ref.C
48mm	43mm	SLE2300

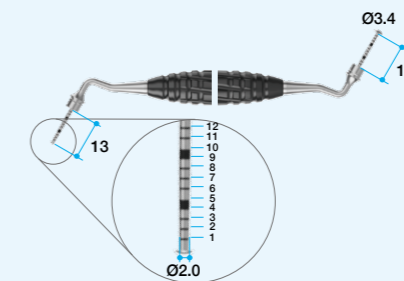

Depth Gauge

Diameter	Ref.C
Ø2.0 / Ø3.4	SCDG2034

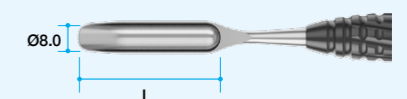

Bone Carrier

Length	Ref.C
35mm	SCBC8035